



with
Japanese
Needle

Elvacryl

POLYGLACTIN 910

Synthetic | Absorbable | Multifilament

KEMENKES RI AKD 21603323442



Elva Surgical Suture

Elvalene
POLYPROPYLENE

Elvadio
POLYDIOXANONE

Elvacryl
POLYGLACTIN 910

Elvamono
POLYGLECAPRONE 25

Elvarapid
POLYGLACTIN 910

Specification	Elvalene	Elvadio	Elvacryl	Elvamono	Elvarapid
Specification	Monofilament, Synthetic	Monofilament, Synthetic	Multifilament, Synthetic	Monofilament, Synthetic	Braided, Synthetic
Color	Blue	Purple	Purple	Violet	Colorless
Complete Absorption	Nonabsorbable	6 months	60 - 70 days	90 - 120 days	56 - 70 days
Tensile Strength (Post Implantation)	Indefinite	2 weeks, approx 75% 4 weeks, approx 65% 6 weeks, approx 50%	2 weeks, approx 75% 3 weeks, approx 50% 4 weeks, approx 25%	gradually	lost by 5 weeks
USP Size	USP 1 - USP 7-0	USP 1 - USP 6-0	USP 2 - USP 5-0	USP 1 - USP 4-0	USP 0 - USP 4-0

Benang Bedah **Elvacryl** adalah benang bedah sintetis, dapat diserap, multifilamen, steril yang terdiri dari kopolimer yang terbuat dari 90% glycolide dan 10% L-lactide. Produk ini dilapisi dengan kopolimer Poly (glycolide-co-lactide) (30/70) dan Calcium Stearate. Produk ini telah terbukti inert, nonantigenic dan nonpyrogenic. Benang Bedah **Elvacryl** tersedia berwarna violet, menggunakan pewarna D&C Violet No. 2.

Benang bedah **Elvacryl** diindikasikan untuk digunakan dalam aproksimasi jaringan kulit dan jaringan lunak secara umum, termasuk penggunaan pada operasi mata. Benang bedah **Elvacryl** tidak digunakan pada jaringan kardiovaskular dan saraf.

Benang bedah **Elvacryl** yang didistribusikan, berwarna dalam ukuran USP 5-0 sampai 2 (metrik 1 sampai 5) dalam panjang yang bervariasi, dengan jarum.

Produk ini dipasarkan dan didistribusikan di bawah entitas **PT Enseval Medika Prima**.

Elvacryl surgical suture is a synthetic, absorbable, multifilament and sterile surgical suture composed of copolymer made from 90% glycolide and 10% L-lactide. Elvacryl is coated with Poly(glycolide-co-lactide) (30/70) copolymer and Calcium Stearate and has been found to be inert, nonantigenic and nonpyrogenic. **Elvacryl** surgical sutures are dyed in violet color, with D&C Violet No.2.

Elvacryl surgical sutures are indicated for use in skin tissue and general soft tissue approximation, including use in ophthalmic surgery. **Elvacryl** surgical sutures are not indicated for use in cardiovascular and neurological tissues.

Elvacryl surgical sutures are supplied in USP sized 5-0 through 2 (metric 1 through 5) in a variety of lengths, with needle.

This product is marketed and distributed under the entity **PT Enseval Medika Prima**.



DESKRIPSI UMUM KODE BENANG BEDAH

GENERAL DESCRIPTION OF SURGICAL SUTURE CODE

Contoh *Example*



1 Inisial Material *Material Initial*

G = Polyglactin 910

2 Ukuran Benang *Thread size*

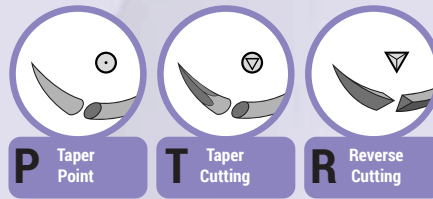
Kode	Ukuran
9	USP 2
1	USP 1
0	USP 0
2	USP 2-0
~	~
5	USP 5-0

4 Panjang Jarum *Needle Length*



Kode = Panjang Jarum
Code = Needle Length

5 Bentuk Ujung Jarum *Point Shape*

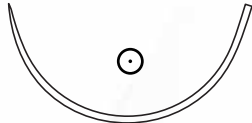

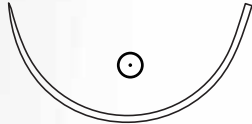



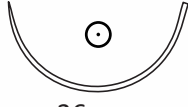


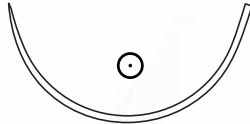



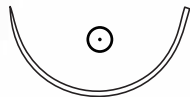


3 Lengkungan *Curvature*








Kode	Lengkungan
1	1/2 180°
2	3/8 135°

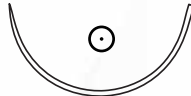


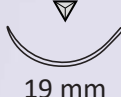
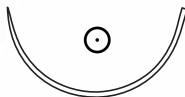


6 Nomor Berjalan *Running Number*








- 00** : Dus isi 12, Panjang Benang 90 cm
Box of contents 12, Suture Length 90 cm
- 01** : Dus isi 18, Panjang Benang 90 cm
Box of contents 12, Suture Length 90 cm
- 02** : Dus isi 24, Panjang Benang 90 cm
Box of contents 12, Suture Length 90 cm



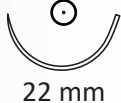

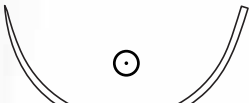


Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G9148P-00	G948P	12	2	 48 mm	1/2 Taper Point
G9148P-01		18			
G9148P-02		24			
G9140P	G940P	12	2	 40 mm	1/2 Taper Point
G9140P-01		18			
G9140P-02		24			
G1148P-00	G148P	12	1	 48 mm	1/2 Taper Point
G1148P-01		18			
G1148P-02		24			
G1140P-00	G140P	12	1	 40 mm	1/2 Taper Point
G1140P-01		18			
G1140P-02		24			
G1140R-00	G140R	12	1	 40 mm	1/2 Reverse Cutting
G1140R-01		18			
G1140R-02		24			
G1136T-00	G136T	12	1	 36 mm	1/2 Taper Cutting
G1136T-01		18			
G1136T-02		24			
G1136P-00	G136P	12	1	 36 mm	1/2 Taper Point
G1136P-01		18			
G1136P-02		24			


Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G0148P-00	G048P	12	0	 48 mm	1/2 Taper Point
G0148P-01		18			
G0148P-02		24			
G0140P-00	G040P	12	0	 40 mm	1/2 Taper Point
G0140P-01		18			
G0140P-02		24			
G0136P-00	G036P	12	0	 36 mm	1/2 Taper Point
G0136P-01		18			
G0136P-02		24			
G2136T-00	G236T	12	2 - 0	 36 mm	1/2 Taper Cutting
G2136T-01		18			
G2136T-02		24			
G2136P-00	G236P	12	2 - 0	 36 mm	1/2 Taper Point
G2136P-01		18			
G2136P-02		24			
G2130R-00	G230R	12	2 - 0	 30 mm	1/2 Reverse Cutting
G2130R-01		18			
G2130R-02		24			
G2126T-00	G226P	12	2 - 0	 26 mm	1/2 Taper Cutting
G2126T-01		18			
G2126T-02		24			

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G2126R-00	G226R	12	2 - 0	 26 mm	1/2 Reverse Cutting
G2126R-01		18			
G2126R-02		24			
G2226R-00	G226R	12	2 - 0	 26 mm	3/8 Reverse Cutting
G2226R-01		18			
G2226R-02		24			
G2126P-00	G226P	12	2 - 0	 26 mm	1/2 Taper Point
G2126P-01		18			
G2126P-02		24			
G2121R-00	G221R	12	2 - 0	 21 mm	1/2 Reverse Cutting
G2121R-01		18			
G2121R-02		24			
G2122T-00	G222T	12	2 - 0	 22 mm	1/2 Taper Cutting
G2122T-01		18			
G2122T-02		24			
G2122R-00	G222R	12	2 - 0	 22 mm	1/2 Reverse Cutting
G2122R-01		18			
G2122R-01		24			
G2122P-00	G222P	12	2 - 0	 22 mm	1/2 Taper Point
G2122P-01		18			
G2122P-02		24			

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G2137P-00	G237P	12	2 - 0	 37 mm	1/2 Taper Point
G2137P-01		18			
G2137P-02		24			
G2140P-00	G240P	12	2 - 0	 40 mm	1/2 Taper Point
G2140P-01		18			
G2140P-02		24			
G3121R-00	G321R	12	3 - 0	 21 mm	1/2 Reverse Cutting
G3121R-01		18			
G3121R-02		24			
G3219R-00	G319R	12	3 - 0	 19 mm	1/2 Reverse Cutting
G3219R-01		18			
G3219R-02		24			
G3136P-00	G336P	12	3 - 0	 36 mm	1/2 Taper Point
G3136P-01		18			
G3136P-02		24			
G3130R-00	G330R	12	3 - 0	 30 mm	1/2 Reverse Cutting
G3130R-01		18			
G3130R-02		24			
G3130P-00	G330P	12	3 - 0	 30 mm	1/2 Taper Point
G3130P-01		18			
G3130P-02		24			

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G3126P-00	G326P	12	3 - 0	 26 mm	1/2 Taper Point
G3126P-01		18			
G3126P-02		24			
G3126T-00	G326T	12	3 - 0	 26 mm	1/2 Taper Cutting
G3126T-01		18			
G3126T-02		24			
G3126R-00	G326R	12	3 - 0	 26 mm	1/2 Reverse Cutting
G3126R-01		18			
G3126R-02		24			
G3125P-00	G325P	12	3 - 0	 25 mm	1/2 Taper Point
G3125P-01		18			
G3125P-02		24			
G3122R-00	G322R	12	3 - 0	 26 mm	1/2 Reverse Cutting
G3122R-01		18			
G3122R-02		24			
G3122T-00	G322T	12	3 - 0	 22 mm	1/2 Taper Cutting
G3122T-01		18			
G3122T-02		24			
G3122P-00	G322P	12	3 - 0	 22 mm	1/2 Taper Point
G3122P-01		18			
G3122P-02		24			

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G4216R-00	G416R	12	4 - 0	 16 mm	3/8 Reverse Cutting
G4216R-01		18			
G4216R-02		24			
G4118P-00	G418P	12	4 - 0	 18 mm	1/2 Taper Point
G4118P-01		18			
G4118P-02		24			
G4122P-00	G422P	12	4 - 0	 22 mm	1/2 Taper Point
G4122P-01		18			
G4122P-02		24			
G5216R-00	G516R	12	5 - 0	 16 mm	3/8 Reverse Cutting
G5216R-01		18			
G5216R-02		24			
G1148P-20	G148P	12	1	 48 mm	1/2 Taper Point
G1148P-21		18			
G1148P-22		24			
G1140P-20	G140P	12	1	 40 mm	1/2 Taper Point
G1140P-21		18			
G1140P-22		24			
G3124R-20	G324R	12	3 - 0	 24 mm	1/2 Reverse Cutting
G3124R-21		18			
G3124R-22		24			

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description
G3224R-00	G324R	12	3 - 0	 24 mm	3/8 Reverse Cutting
G3224R-01		18			
G3224R-02		24			



Manufactured by :

PT Forsta Kalmedic Global

A Kalbe Company

Jl. Pulogadung No.23, Industrial Estate (JIEP), Jakarta 13930

Phone : 021 - 5086 7788

<https://forsta.co.id>



Distributed by :

PT Enseval Medika Prima

A Kalbe Company

Jl. Pulo Lentut No.12, Industrial Estate (JIEP), Jakarta 13930

Phone : 021 - 4682 3234

<https://emp.co.id>