



with Japanese Needle

Elvamono

POLYGLECAPRONE 25

Synthetic | Absorbable | Monofilament

KEMENKES RI AKD 21603323371

 **Elva Surgical Suture**
INDONESIA'S FIRST PREMIUM SUTURE



Elva Surgical Suture

Elvalene
POLYPROPYLENE

Elvadio
POLYDIOXANONE

Elvacryl
POLYGLACTIN 910

Elvamono
POLYGLECAPRONE 25

Elvarapid
POLYGLACTIN 910

Specification	Elvalene	Elvadio	Elvacryl	Elvamono	Elvarapid
Specification	Monofilament, Synthetic	Monofilament, Synthetic	Multifilament, Synthetic	Monofilament, Synthetic	Braided, Synthetic
Color	Blue	Violet	Violet	Violet	Colorless
Complete Absorption	Nonabsorbable	6 months	60 - 70 days	90 - 120 days	56 - 70 days
Tensile Strength (Post Implantation)	Indefinite	2 weeks, approx 75% 4 weeks, approx 65% 6 weeks, approx 50%	2 weeks, approx 75% 3 weeks, approx 50% 4 weeks, approx 25%	gradually	lost by 5 weeks
USP Size	USP 1 - USP 7-0	USP 1 - USP 6-0	USP 2 - USP 5-0	USP 1 - USP 4-0	USP 0 - USP 4-0

Benang Bedah **Elvamono** adalah benang bedah sintesis yang dapat diserap, monofilamen, terdiri dari dari kopolimer yang terbuat dari glikolida dan epsilon-kaprolakton. Produk ini telah terbukti non-antigenik, non-pirogenik dan hanya menimbulkan sedikit reaksi selama penyerapan. Benang Bedah **Elvamono** tersedia berwarna Violet.

Benang Bedah **Elvamono** diindikasikan untuk digunakan dalam aproksimasi jaringan lunak secara umum dan/atau ligasi. Benang Bedah **Elvamono** tidak digunakan dalam jaringan kardiovaskular atau neurologis, bedah mikro atau bedah mata.

Benang Bedah **Elvamono** didistribusikan dalam kondisi steril dalam ukuran USP dan variasi panjang benang, beserta jarum. Tersedia dalam kemasan isi 12, 18 dan 24 Benang Bedah **Elvamono**.

Produk ini dipasarkan dan didistribusikan di bawah entitas **PT Enseval Medika Prima**.

Elvamono Surgical Suture is an absorbable, monofilament, synthetic surgical suture composed of a copolymer made of glycolides and epsilon-caprolactone. This product is proven to be non-antigenic, non-pyrogenic and causes little reaction during absorption. **Elvamono** Surgical Suture is available in Violet.

Elvamono Surgical Suture is indicated for use in general soft tissue approximation and/or ligation. **Elvamono** Surgical Thread is not used in cardiovascular or neurological tissue, microsurgery or eye surgery.

Elvamono Surgical Sutures are distributed sterile with varying USP sizes and thread lengths, along with needles. Available in packaging of 12, 18, and 24 **Elvamono** Surgical Sutures.


This product is marketed and distributed under the entity **PT Enseval Medika Prima**.



DESKRIPSI UMUM KODE BENANG BEDAH

GENERAL DESCRIPTION OF SURGICAL SUTURE CODE

Contoh Example

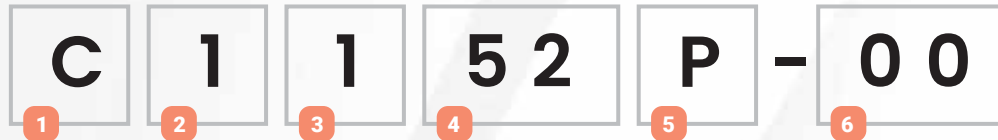


1 Inisial Material *Material Initial*

C = Polyglecaprone 25

2 Ukuran Benang *Thread size*

Kode	Ukuran
1	USP 1
0	USP 0
2	USP 2-0
3	USP 3-0
4	USP 4-0

3 Lengkungan *Curvature*


Kode	Lengkungan
1	1/2 
2	3/8 

4 Panjang Jarum *Needle Length*



Kode = Panjang Jarum
Code = Needle Length

5 Bentuk Ujung Jarum *Point Shape*



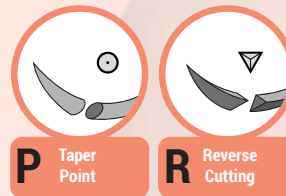








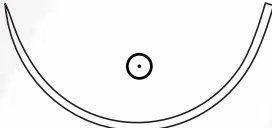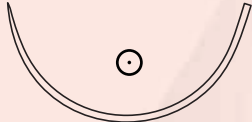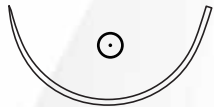

P Taper Point

R Reverse Cutting

6 Nomor Berjalan *Running Number*

- 00** : Dus isi 12, Panjang Benang 90 cm
Box of contents 12, Suture Length 90 cm
- 01** : Dus isi 18, Panjang Benang 90 cm
Box of contents 12, Suture Length 90 cm
- 02** : Dus isi 24, Panjang Benang 90 cm
Box of contents 12, Suture Length 90 cm

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description	Indication
C1152P-00	C152P	12	1		1/2 Taper Point	Muscle, fascia, uterus, general closure
C1152P-01		18				
C1152P-02		24				
C1148P-00	C148P	12	1		1/2 Taper Point	Muscle, fascia, obgyn, general closure
C1148P-01		18				
C1148P-02		24				
C1140P-00	C140P	12	1		1/2 Taper Point	Muscle, fascia, uterus, general closure
C1140P-01		18				
C1140P-02		24				
C0140P-00	C040P	12	0		1/2 Taper Point	Muscle, fascia, obgyn, general closure
C0140P-01		18				
C0140P-02		24				
C0148P-00	C048P	12	0		1/2 Taper Point	Muscle, fascia, obgyn, general closure
C0148P-01		18				
C0148P-02		24				
C0137P-00	C037P	12	0		1/2 Taper Point	Muscle, fascia, obgyn, general closure
C0137P-01		18				
C0137P-02		24				
C2148P-00	C248P	12	2 - 0		1/2 Taper Point	Muscle, fascia, obgyn, general closure
C2148P-01		18				
C2148P-02		24				

Principal Code	EMP Code	Box of Content	USP Size	Needle Length	Needle Description	Indication
C3219R-00	C319R	12	3 - 0	 19 mm	3/8 Reverse Cutting	Skin, ubcutaneous, subcuticular
C3219R-01		18				
C3219R-02		24				
C3224R-00	C324R	12	3 - 0	 24 mm	3/8 Reverse Cutting	Skin, ubcutaneous, Plastic, subcuticular
C3224R-01		18				
C3224R-02		24				
C3130R-00	C330R	12	3 - 0	 30 mm	1/2 Reverse Cutting	Skin, ubcutaneous, subcuticular
C3130R-01		18				
C3130R-02		24				
C4219R-00	C419R	12	4 - 0	 19 mm	3/8 Reverse Cutting	Skin, ubcutaneous, subcuticular
C4219R-01		18				
C4219R-02		24				



Manufactured by :

PT Forsta Kalmedic Global

A Kalbe Company

Jl. Pulogadung No.23, Industrial Estate (JIEP), Jakarta 13930

Phone : 021 - 5086 7788

<https://forsta.co.id>



Distributed by :

PT Enseval Medika Prima

A Kalbe Company

Jl. Pulo Lentut No.12, Industrial Estate (JIEP), Jakarta 13930

Phone : 021 - 4682 3234

<https://emp.co.id>